

PRO CART II™

Mobile Treatment Console




DNTLworks™
Equipment Corporation

THE INNOVATIVE PRO CART II™ COMBINES THE BEST ATTRIBUTES OF SELF-CONTAINED MOBILITY AND A HIGH-TECH DESIGN THAT IS PERFECTLY SUITED FOR PROFESSIONAL IN-OFFICE USE.

FEATURES INCLUDE:

- ♥ *Sophisticated, one-piece design: completely self-contained with built-in air compressor, separate vacuum pump, water & waste management systems.*
- ♥ *High-strength aluminum chassis with innovative WhisperQuiet™ sound reduction technology.*
- ♥ *Separate oil-free adjustable vacuum pump for powerful suction.*
- ♥ *Large, four liter waste container with automatic overflow shutoff and electric waste purge.*
- ♥ *Ample, 1.5 liter water reservoir with additional external water connection.*
- ♥ *ARL* Certified for compliance with Safety Standard UL 544 "Medical & Dental Equipment"*

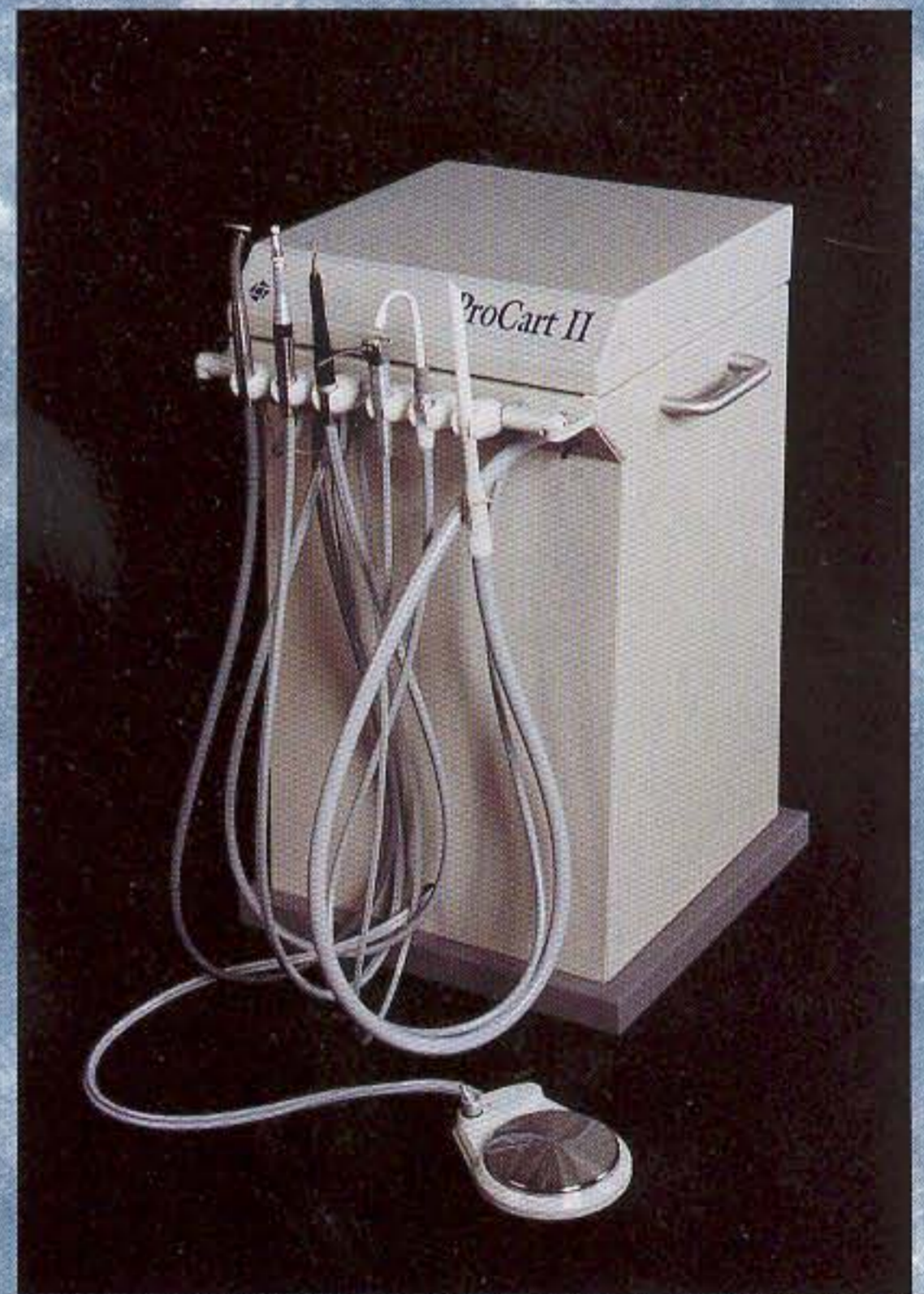

Designed specifically for use in the modern dental office, nursing home, hospital and medical center, the ProCart II combines superior performance with state-of-the-art design.

ProCart II, WhisperQuiet and DNTLworks are trademarks of
DNTLworks Equipment Corporation © 1998



*Applied Research Laboratories (ARL) is a nationally recognized independent testing, listing, labeling and follow-up service laboratory.



PRO CART II™

Mobile Treatment Console



DNTLworks™

SPECIFICATIONS

DIMENSIONS/WEIGHT

Width: 20"
Depth: 21"
Height: 30" to 34" adjustable
Weight: 100 lbs

ELECTRICAL

Power Supply: 115 Vac, 60 Hz or
220 Vac, 50 Hz
Circuit breaker protection: 15 amps
Hospital grade wiring, NEMA type
5-15P Plug

COMPRESSOR

Oil-free 1/3 hp motor
Working pressure: 70 to 100 psi
Flow rate: 2.5 cfm
Sound level: 48 dB
Air reservoir: 4 liters
Automatic air tank purge

VACUUMMOTOR

Oil-free 1/3 hp motor
Adjustable vacuum: 0" to 21" Hg
Vacuum flow: 4.6 cfm
Waste receptacle: 4 liters
Sound level: 48 dB
Automatic overflow shutoff
Aseptic filtered exhaust
Large volume solids trap
Automatic waste removal

WATERMODULE

1.5 liter reservoir
Pressurized system: 30 to 40 lbs
Aseptic non-retracting HPC coolant
Handpiece flush system
External quick connects for
water and air



***Sleek Design for
In-office Use***



***Easy Access to
Control Panel***



***Auxillary Air & Water
Quick Connects;
Automatic Waste Removal;
Additional Electrical Outlet***



***Swing-Out Instrument
Arms***

- *Wet/dry foot control*
- *Automatic controls for one, two and three handpieces*
- *Accommodates most brands of handpieces and fiber optic systems*
- *Three-way syringe*
- *Adjustable H.V.E. system*
- *Adjustable saliva ejector*
- *Asepsis designed chassis and handpiece holders*
- *Recessed controls for easy cleaning*
- *Air filter moisture separator*
- *Repair/Tool kit with Operations Manual*
- *Side maneuvering handles and high-grade casters for easy movement*
- *Optional arm-mounted trays available*
- *Optional cover available*

Lions Dental Supply
www.LionsDentalSupply.com
951-276-3225