



LITEX™ 695C

TWO-PHASE POLYMERIZATION CORDLESS LED CURING LIGHT

LITEX™
695C

- Thin and lightweight dual intensity light polymerization LED curing light.
- Patent pending 360° rotatable light guide with front loaded LED provides easy access and maximum power.
- No more bulbs to replace with average LED life of 50,000 hours.
- Smooth and contoured ergonomic design offers added comfort and easy cleaning.
- Conveniently located thumb buttons set timer and intensity for two stage step-curing.
- Two stage polymerization with initial lower intensity prior to higher intensity polymerization resulting in less shrinkage stress.
- Automatic timer with audible tone and digital display on handpiece.
- Light intensity of 1,200 mW/cm² cures composite in 10 seconds / 4 mm.
- Curing cycle can be interrupted manually with one-touch trigger button.
- Large lithium battery capacity provides 750 ten-second curing cycles on a full charge.
- Low battery warning beeper.
- Built-in light intensity indicator.
- Eye protection shield included for safe operation.

ISO 9001
ISO 13485




GMP

DENTAMERICA®

LITEX™ 695C

TWO-PHASE POLYMERIZATION CORDLESS LED CURING LIGHT

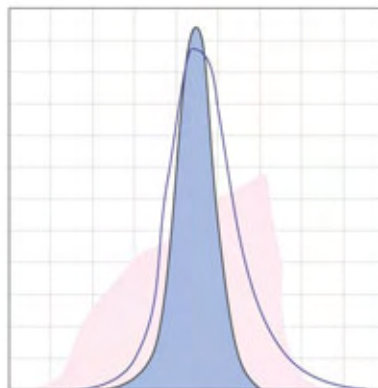

nm 400 420 440 460 480 500 520 540

■ LITEX 695C LED
— Camphorquinone
■ Halogen Lamps

10-second cure for 4mm curing depth.



Built-in Light Intensity Indicator.



Reorder No.69502
Blue LED Probe



Reorder No.69504
Diagnostic White LED Probe

Patent pending 360°rotatable light guide with front loaded LED provides easy access and maximum power



Reorder No.664
Eye Protective Shield



Reorder No.69506
LED Probe Cover
For LED protection

SPECIFICATIONS:

- INPUT: 100-240V AC 0.4A
- OUTPUT: 9V DC 1.3A
- HANDPIECE: 130g
- CHARGING UNIT: 260g
- L170xW38.5xH33mm
- L147xW93xH150mm

DENTAMERICA®

DISTRIBUTOR

Lions Dental Supply
951-276-3225
www.LionsDentalSupply.com

DENTAMERICA®