


Blendex™


Automatic Alginate Mixer



① Add alginate and water to the mixing cup.

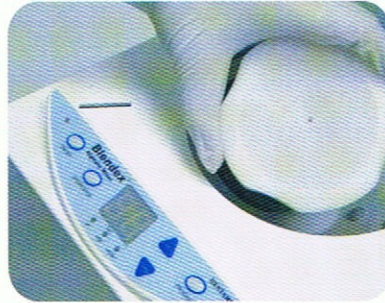


② Place the mixing cup inside well.

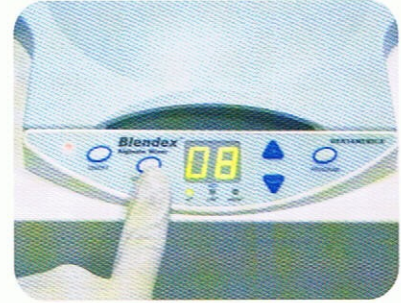

③ Push the START button for a perfectly smooth **Blendex™** mix.

3 customizable memory presets available (settings can be changed)

Cups	Mixing Time (sec)
	8
	10
	12

SPECIFICATIONS

Dimensions	W205 x D245 x H295 (mm)
Weight	18 kg
Voltage	AC220V~230V (50/60 Hz) AC100V~120V (50/60 Hz)
Maximum Power	400W
Timer	1-16 sec, adjustable by time
Memory	3 (installed settings: 8, 10, 12 sec)
Frequency	3,600 RPM

